

NR-fit Anaesthesia Product Catalogue



Shanghai Kindly Enterprise Development Group Co., Ltd.

Add: No.658 Gaochao Road, Shanghai,China

Post Code: 201803

Tel: 86-21-69119653 59144298

Fax: 86-21-59118500

Http: // www.kdlchina.com

Email: webmaster@kdlchina.com

Rev . : 20191030

SHANGHAI KINDLY
ENTERPRISE DEVELOPMENT GROUP CO.,LTD.

To address the risk of wrong route administration ISO 80369-6 medical device connectors are being progressively introduced worldwide, replacing conventional Luer connectors. Sampling and delivery devices intended for neuraxial applications and for major regional anaesthesia will in future utilise ISO 80369-6 connectors.

Neuraxial connectors devices will for easy identification be coloured yellow.

Product Catalogue

NR-fit Syringe -----01

Anaesthesia Needles -----02

NR-fit Stopcocks -----03

NR-fit Extension Sets -----04



NR-fit Syringe

Structure and composition:

The NR-fit syringes are assembled by Barrel, Plunger, Piston. All parts and material for this product meet medical requirements, after sterilized by ETO, pyrogen free.

Material:

PP, Isoprene rubber, Silicone oil



Remark: The type can make as clients' request.


Anaesthesia Needles

Structure and composition:

The Anaesthesia needles consist of needle hub, needle tube (out), needle tube (inner), pro tector cap. Major categories: Spinal anesthesia needles, Epidural anesthesia needles, Combined anaesthesia needles.

Material:

PP,PC,ABS, SUS304 Stainless Steel Cannula,Silicone oil.

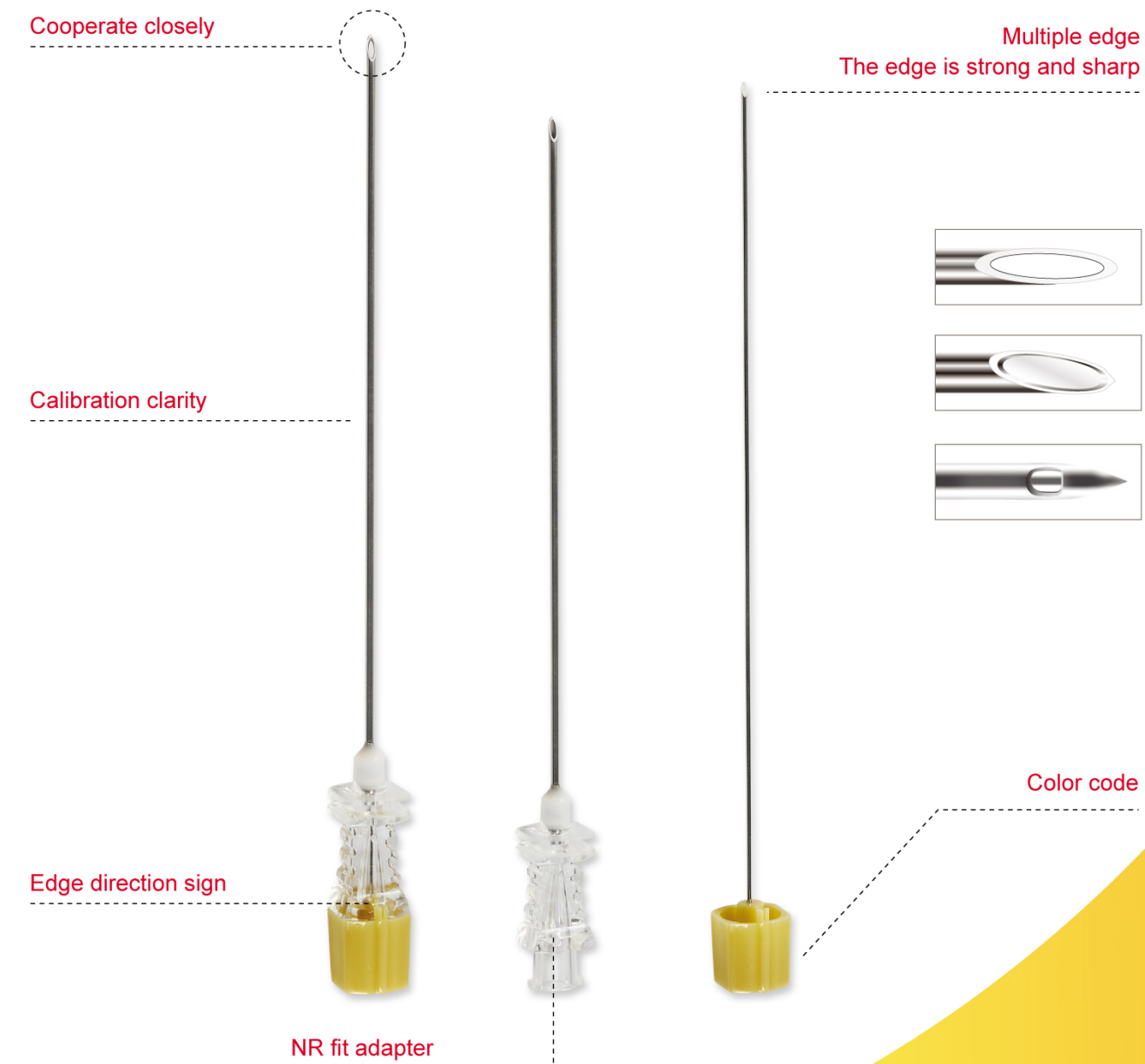


NR-fit Stopcocks

Structure and composition:

The NR-fit Stopcocks are assembled by stopcock body, Locking screw. All parts of this product and processing of materials are used in line with the requirements of medical materials, after sterilization by EO, pyrogen free.

Material:

PCTG, HDPE.

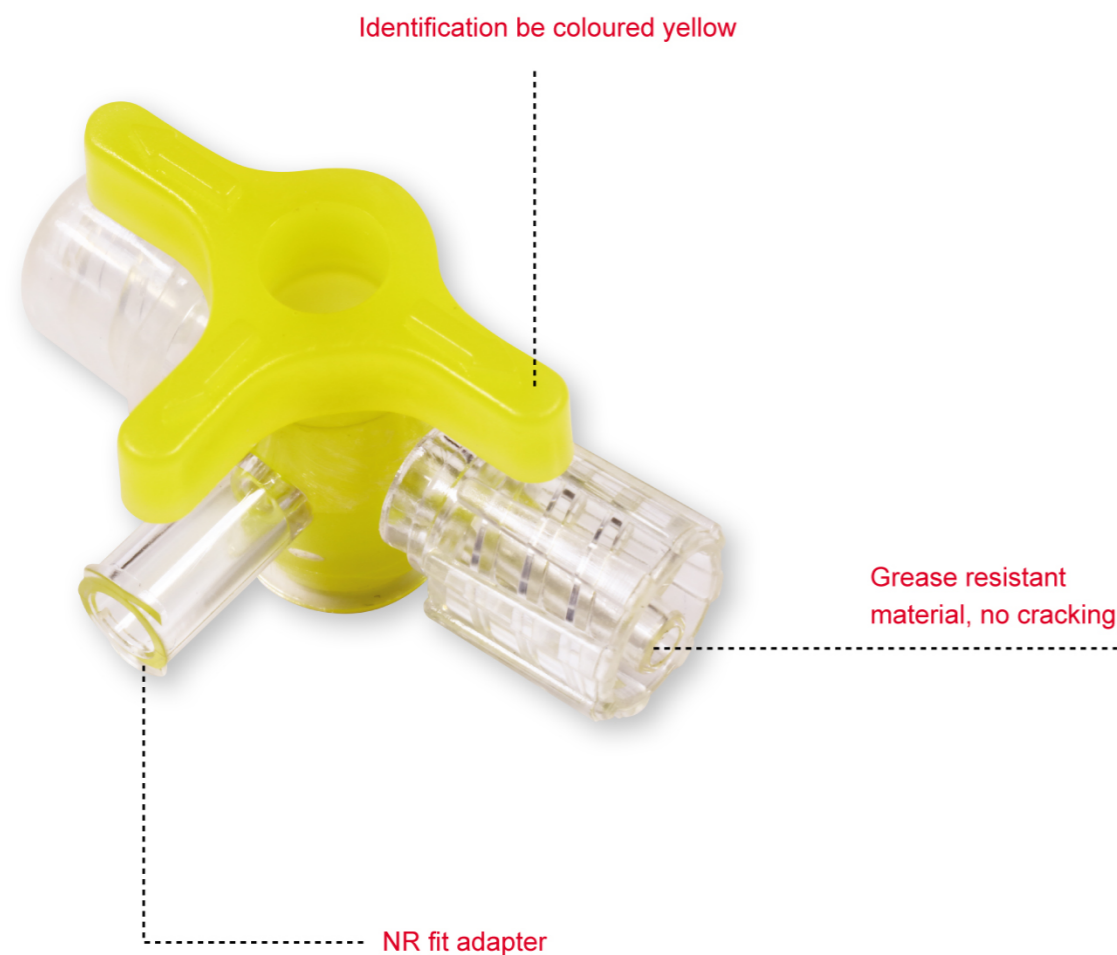

Remark: The type can make as clients' request.

NR-fit Extension Sets

Structure and composition:

The NR-fit Extension sets are assembled by NR-fit adapter, Robert lock clamp, tube. All parts of this product and processing of materials are used in line with the requirements of medical materials, after sterilization by EO, pyrogen free.

Material:

PCTG, PVC, POM/HDPE.

